hâws  
RELAX

30-HSRB21

## RECOVERY BOOT

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Bruksanvisning - Svenska

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Instruction manual - English

CE

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## EN - INSTRUCTION MANUAL RECOVERY BOOT

Please read this manual carefully before using the product and keep it in a safe place for future reference.

**Note:** The actual appearance of the product may differ slightly from illustrations.

### PRECAUTIONS

1. Read the instructions carefully before use. Do not use the product for extended periods to avoid muscle discomfort or strain.
2. Only qualified service personnel may disassemble or repair this product. Unauthorised repairs are prohibited.
3. If the device is damaged, discontinue use and have it repaired by the manufacturer or an authorised service centre to avoid risk.
4. Do not use the product in humid environments or if it has come into contact with water or other liquids.
5. After use, store the device in a dry, cool place away from direct sunlight.
6. Keep away from heat sources, flammable substances, chemical solutions, or high-humidity areas.
7. Avoid contact with sharp objects to prevent damage to the device.
8. If you experience any discomfort or unusual symptoms during use, stop using the device immediately.
9. Adjust usage duration and intensity according to personal comfort to avoid overuse or strain.

#### **Consult a doctor or use under professional supervision if you:**

- Are pregnant (especially high-risk pregnancies)
- Are a child
- Have infectious skin conditions, redness, swelling, or burns
- Have a serious medical condition that contraindicates massage
- Have weak or damaged bones
- Use a pacemaker
- Have metal implants, such as prosthetics
- Suffer from severe heart disease or abnormal blood pressure

**Product Protection and Storage:**

1. Keep the device out of reach of children.
2. Do not place any objects on the main unit or the fabric cover.
3. This product is designed only for use on the legs and feet. Do not use it on other parts of the body.
4. Use only the charging cable provided with the product and a compatible, certified adapter.
5. The device must not be used while charging.
6. The fabric areas are not washable. To clean, use a damp, well-wrung cloth. Use a soft cloth to clean the main unit.
7. Before first use—or if the battery has been fully discharged after long-term storage—please charge the product fully.
8. If not in use for an extended period, switch off the device and recharge it every two months to preserve battery life.


## PRODUCT DESCRIPTION

### Technical Specifications

Charging voltage:	5 V
Power consumption:	10 W
Battery capacity:	2000 mAh

### Overview

Main unit
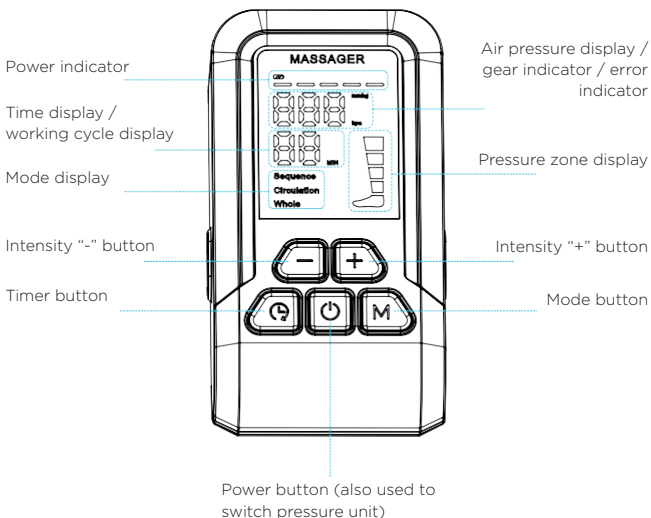
## PRODUCT PRESENTATION

- Power Button**
  - Press and hold to switch the device on or off.
  - Press to toggle the air pressure unit between kPa and mmHg.
- Timer Button**
  - Press to cycle through timer settings: 10, 20, or 30 minutes.
- Mode Button**
  - Press to switch between modes 1, 2, and 3.
  - (Note: The massage function pauses briefly when switching modes.)
  - Press and hold to enter or exit fixed-point massage mode.
- Intensity Buttons**
  - Press the “-” button to reduce massage intensity.
  - Press the “+” button to increase massage intensity.

### Massage Modes

There are three intelligent massage modes:

#### Mode 1 - Sequential

Designed to relax muscles with individual inflation of air bags in a sequential pattern.

#### Mode 2 - Circulation

Aims to restore physical energy using simultaneous inflation of two air bags.

#### Mode 3 - Overall

Provides full care by combining inflation of one, two, or all three air bags in various patterns.

## Fixed-Point Massage Selection

- To enter fixed-point massage mode, press and hold the mode button until you hear two beeps.
- Once active, press the mode button to cycle through the air bags you wish to target.
- The device will automatically start the fixed-point massage 5 seconds after selection.
- To exit this mode, press and hold the mode button again.

## Adjustable Pressure Levels

Each mode supports five levels of pressure for a personalised massage experience.

## Working Cycle

1. **Timing Options:**
  - The default runtime is 20 minutes upon startup.
  - Press the timer button to cycle through 30 minutes → 20 minutes → 10 minutes.
2. **Quick Removal:**
  - In case of emergency or if you wish to stop use, you can quickly remove the main unit from the fabric sleeve.

## Power Display

- The power level is indicated by five green LED lights.
- When the battery is low, the battery icon will flash.
- Charging status is also shown via the power indicator.  
**Note:** The device cannot be used while charging.

## Charging Instructions

- Switch off the main unit before charging and remove it from the sleeve.
- Use the USB-C charging port located on the main unit.
- When removing the device, press the indented grips on both sides simultaneously to release it.
- When reinstalling, align the main unit with the installation slot and ensure it is securely fitted.

## PRECAUTIONS & MAINTENANCE

If you experience any issues while using the product, please refer to the following troubleshooting tips. If the problem persists or the device is not functioning properly, switch off the power immediately.

**Do not attempt to disassemble or repair the product yourself.** Unauthorised repairs will void the product warranty. The manufacturer accepts no responsibility for damage resulting from improper handling.

### Troubleshooting Guide

Problem	Possible Cause / Solution
Product does not power on	Check the battery and ensure the power button is pressed correctly.
	Ensure the massage button has been pressed and that the massage head is securely attached.
Product is not charging or displaying	Check that the USB-C charging cable is properly connected to the charging port.
Product stops during operation	Check if the battery is low or if the automatic timer has reached the set time limit.
	If the problem continues, contact the authorised service centre.

### Danger

- This product uses a specialised rechargeable battery that must not be used in other devices.
- Do not attempt to recharge a battery once it has been removed from the device.
- Do not throw the battery into fire or expose it to high heat.
- Do not puncture, crush, or dismantle the battery.
- Do not allow metal objects to come into contact with the battery terminals.
- Keep away from metal items such as necklaces and hairpins.
- Do not charge, use, or store the battery near heat sources or in direct sunlight.
- Exposure to high temperatures may cause overheating, fire, or explosion.

## Warning

- Keep removed batteries out of the reach of infants and children. If swallowed, seek immediate medical attention.
- If battery fluid leaks:
  - **Eyes:** Do not rub. Rinse thoroughly with clean water and seek medical help immediately – risk of blindness.
  - **Skin or clothing:** Rinse with clean water. If irritation persists, seek medical advice.
- Improper handling may result in fire, electric shock, or injury.

### Protecting the environment

This symbol indicates the separate collection of electrical and electronic equipment. The main objective of the separate collection of this waste is to reduce the quantity to be disposed of, promote re-use, recycling and other forms of recovery in order to reduce its negative effects on the environment. When it is no longer in use, this appliance must not be disposed of together with other unsorted municipal waste. The user is responsible for delivering it free of charge to the selective collection facilities available for that purpose. Pieces of equipment that are no longer in use can also be taken back at points of sale, by purchasing new pieces of equipment that are equivalent and have the same functions. For more detailed information on collection points, please contact your City Council or a point of sale for this equipment.

**FUN NORDIC ApS**

Rugkobbel 260

6200 Aabenraa

Denmark

[www.funnordic.dk](http://www.funnordic.dk)

[info@funnordic.dk](mailto:info@funnordic.dk)